

# TREATING THE COMPLETE DENTURE PATIENT

EDITED BY

CARL F. DRISCOLL | WILLIAM GLEN GOLDEN



WILEY Blackwell



## Treating the Complete Denture Patient



## Treating the Complete Denture Patient

*Edited by*

*William Glen Golden, DDS*

*Associate Clinical Professor (retired)*

*The Ohio State University College of Dentistry, Columbus, OH, USA*

*assisted by*

*Carl F. Driscoll, DMD*

*Professor*

*School of Dentistry, University of Maryland, Baltimore, MD, USA;*

*Private Practice, Bethesda, MD, USA*

**WILEY** Blackwell

This edition first published 2020  
© 2020 John Wiley & Sons, Inc.

All rights reserved. No part of this publication may be reproduced, stored in a retrieval system, or transmitted, in any form or by any means, electronic, mechanical, photocopying, recording or otherwise, except as permitted by law. Advice on how to obtain permission to reuse material from this title is available at <http://www.wiley.com/go/permissions>.

The rights of Carl F. Driscoll and William Glen Golden to be identified as the authors of the editorial material in this work has been asserted in accordance with law.

*Registered Office*

John Wiley & Sons, Inc., 111 River Street, Hoboken, NJ 07030, USA

*Editorial Office*

111 River Street, Hoboken, NJ 07030, USA

For details of our global editorial offices, customer services, and more information about Wiley products visit us at [www.wiley.com](http://www.wiley.com).

Wiley also publishes its books in a variety of electronic formats and by print-on-demand. Some content that appears in standard print versions of this book may not be available in other formats.

*Limit of Liability/Disclaimer of Warranty*

The contents of this work are intended to further general scientific research, understanding, and discussion only and are not intended and should not be relied upon as recommending or promoting scientific method, diagnosis, or treatment by physicians for any particular patient. In view of ongoing research, equipment modifications, changes in governmental regulations, and the constant flow of information relating to the use of medicines, equipment, and devices, the reader is urged to review and evaluate the information provided in the package insert or instructions for each medicine, equipment, or device for, among other things, any changes in the instructions or indication of usage and for added warnings and precautions. While the publisher and authors have used their best efforts in preparing this work, they make no representations or warranties with respect to the accuracy or completeness of the contents of this work and specifically disclaim all warranties, including without limitation any implied warranties of merchantability or fitness for a particular purpose. No warranty may be created or extended by sales representatives, written sales materials or promotional statements for this work. The fact that an organization, website, or product is referred to in this work as a citation and/or potential source of further information does not mean that the publisher and authors endorse the information or services the organization, website, or product may provide or recommendations it may make. This work is sold with the understanding that the publisher is not engaged in rendering professional services. The advice and strategies contained herein may not be suitable for your situation. You should consult with a specialist where appropriate. Further, readers should be aware that websites listed in this work may have changed or disappeared between when this work was written and when it is read. Neither the publisher nor authors shall be liable for any loss of profit or any other commercial damages, including but not limited to special, incidental, consequential, or other damages.

*Library of Congress Cataloging-in-Publication Data*

Names: Driscoll, Carl F., editor. | Golden, William Glen, 1947- editor.

Title: Treating the complete denture patient / edited by Carl F. Driscoll, William Glen Golden.

Description: Hoboken, NJ : Wiley-Blackwell, 2020. | Includes bibliographical references and index.

Identifiers: LCCN 2019039403 (print) | LCCN 2019039404 (ebook) | ISBN 9781119569589 (hardback) |

ISBN 9781119569572 (adobe pdf) | ISBN 9781119569565 (epub)

Subjects: MESH: Denture, Complete

Classification: LCC RK656 (print) | LCC RK656 (ebook) | NLM WU 530 | DDC 617.6/92-dc23

LC record available at <https://lcn.loc.gov/2019039403>

LC ebook record available at <https://lcn.loc.gov/2019039404>

Cover Design: Wiley

Cover Image: © hoozone/Getty Images

Set in 9.5/12.5pt STIXTwoText by SPi Global, Pondicherry, India

## Contents

Foreword *vii*

Acknowledgements *ix*

About the Companion Website *xi*

- 1 The Diagnostic Appointment *1*
- 2 Preliminary Impressions *5*
- 3 Custom Trays: The Key to a Great Final Impression *17*
- 4 Using Custom Trays to Make Final Impressions *25*
- 5 The Final Impression *41*
- 6 Boxing and Pouring a Final Impression *47*
- 7 Base Plates and Occlusion Rims *51*
- 8 Adjusting Wax Rims to Fit the Patient *57*
- 9 Using the Denar Facebow *59*
- 10 Mounting Maxillary Casts on an Articulator *63*
- 11 Centric Relation Records *67*
- 12 Selecting Proper Denture Teeth *73*
- 13 Setting Anterior Denture Teeth *81*
- 14 Setting Posterior Anatomic Teeth in Bilaterally Balanced Occlusion *85*
- 15 Setting Zero-Degree Posterior Teeth in Monoplane Occlusion *93*
- 16 Overview of Lingualized Occlusion *99*
- 17 Setting Teeth in a Lingualized Occlusion *105*
- 18 The Anterior Esthetic Wax Try-In *111*
- 19 The Complete Denture Wax-Up *115*

20	Sealing Complete Dentures to the Cast	121
21	The Wax Try-in of Teeth	123
22	Making an Arbitrary Posterior Palatal Seal in the Lab	135
23	Interacting with Dental Labs	137
24	Flasking, Packing, and Processing Complete Dentures	139
25	Completing a Laboratory Remount	151
26	Making a Remount Index	155
27	Remove Complete Denture from the Cast and Finish	159
28	Quality Standards in Appearance	165
29	Inserting Complete Dentures in the Mouth	169
30	Performing a Clinical Remount	175
31	Equilibrating Dentures Set in Lingualized Occlusion	187
32	Equilibrating Fully Anatomic (30°) Set-Ups	191
33	Equilibrating Complete Dentures Set in a Monoplane Occlusion	197
34	Troubleshooting Complete Denture Problems	199
35	Replacing Teeth on a Complete Denture	215
36	Repairing a Broken Complete Denture	219
37	Using Implants to Stabilize a Complete Denture	229
38	Immediate Complete Dentures	239
39	Tissue Conditioners	253
40	The Fabrication of Digital Complete Dentures	263

List of Captions 271

Index 293

## Foreword

**William Glen Golden** began his dental training with nine years enlisted service as a dental prosthetic laboratory technician in the US Navy, achieving the rank of Dental Technician Chief.

He completed dental technician A&C schools in San Diego and Field Medical Service School at Camp Pendleton, California. He served as a dental prosthetic technician and dental clinic supervisor before entering dental school in 1976 at West Virginia University on the Navy version of the Armed Forces Health Professions Scholarship. Graduating in 1980, he was sent to the Naval Dental Clinic in Yokosuka, Japan, where he served for four years, part of that time on the USS *Blue Ridge* (LCC-19), the flagship of the Seventh Fleet.

He completed a postdoctoral Fellowship in Prosthodontics at the Naval Dental Clinic, Great Lakes, Illinois, in 1985 and a Residency in Prosthodontics at the US Navy Dental Postgraduate School in Bethesda, Maryland, in 1988 with a certificate in Prosthodontics.

Following graduation, he served aboard the USS *Nimitz* (CVN 68), at the US Naval Academy in Annapolis, Maryland, as Head of the Prosthodontics Department, as a prosthodontist at the Washington Navy Yard, and as the Director of the Area Dental Prosthetic Laboratory at Naval Dental Clinic, Great Lakes, until his retirement in 1995.

After retiring from the US Navy, Dr Golden worked for two years as a restorative dentist and prosthodontist for a managed-care dental practice in Green Bay, Wisconsin, then as Clinical Assistant Professor and Prosthodontist at the Ohio State University (OSU) College of Dentistry in November 1997. He directed three undergraduate complete denture courses and was promoted to Clinical Associate Professor. He was serving as the Director of Removable Prosthodontics and Director of the Complete Denture Clinic at OSU when he retired at the end of December in 2010.

Following his retirement from OSU, Dr Golden accepted a part-time position as a Clinical Associate Professor at West Virginia University School of Dentistry in November 2011. He became Course Director of a six credit-hour postdoctoral complete denture course and served for two

years until he retired from actively practicing dentistry in June 2014.

**Carl F. Driscoll** is an acclaimed educator, researcher, and clinician. He currently serves as a Professor at the University of Maryland Dental School, where he previously served as Director of the Prosthodontic Residency for 21 years. That assignment followed a 20-year career with the US Army culminating with being the Prosthodontic Program Director at Walter Reed Army Medical Center for three years. Dr Driscoll has served as President of the American Academy of Fixed Prosthodontics, President of the American Board of Prosthodontics, and President of the American College of Prosthodontists. He has been awarded both the Garver-Staffanou Award (for outstanding service as a program director) and the Moulton Award (for outstanding contributions to fixed prosthodontics) from the American Academy of Fixed Prosthodontics. In addition, he recently was awarded the Educator of the Year Award from the American College of Prosthodontists. Dr Driscoll has given over 450 presentations nationally and internationally and has authored over 75 publications. Besides teaching in the PG Pros program at Maryland, he also maintains a private practice in Bethesda, Maryland, with his wife, Dr Sarit Kaplan.

**Nadim Z. Baba** received his DMD degree from the University of Montreal in 1996. He completed a Certificate in Advanced Graduate Studies in Prosthodontics and a Masters degree in Restorative Sciences in Prosthodontics from Boston University School of Dentistry in 1999. Dr Baba is a Professor in the Advanced Education program in Prosthodontics at Loma Linda University School of Dentistry, an Adjunct Professor at the University of Texas Health Science Center School of Dentistry in the Comprehensive Dentistry Department, and maintains a part-time private practice in Glendale, CA. He is currently the President of the American College of Prosthodontists and an active member of various professional organizations, a Diplomate of the American Board of Prosthodontics, and a Fellow of the American College of Prosthodontists and the Academy of Prosthodontics.



## Acknowledgements

To Susan G. Kestner, who guided me in the right direction in organizing my thoughts and putting them down in print.

To Dr. Charles Goodacre, Dr. Timothy Saunders, and Dr. Alejandro Peregrina who encouraged me to write this book, reviewed its content, and made suggestions in how to improve it.

To Anthony Buffamonte, whose lab expertise, encouragement, and friendship I have treasured.

To Dr. Stephen Ancowitz, whose praise and support were always greatly appreciated.

To CAPT Robert Slater, DC, USN, CAPT Charles R. Linkenbach, DC, USN, RADM George A. Besbekos, DC,

USN, and CAPT Alexander Sanderson, DC, USN who recognized my competence as a Dental Prosthetic Technician and greatly encouraged me to become a dentist, then supported me greatly to be selected for the Armed Forces Health Professions Scholarship Program.

To Dr. Carl F. Driscoll, who agreed to be co-author and advisor and help me get this book published.

To Dr. Nadim Z. Baba who agreed to write the chapter on digital complete dentures.

And finally, to my wife, Louise Teets Golden, who put up with the long hours it took me to complete this textbook and kept my schedule organized.



## About the Companion Website

Don't forget to visit the companion website for this book.

[www.wiley.com/go/driscoll/denture](http://www.wiley.com/go/driscoll/denture)



There you will find valuable video material designed to enhance your learning.  
Scan this QR code to visit the companion website





# 1

## The Diagnostic Appointment

Our goal is to teach students what to look for when diagnosing conditions that impact upon the use and prognosis of a complete denture, so that they may be better able to provide a complete denture service to their patients. First, certain anatomical features will be evaluated, and a judgment rendered as to their possible effects on the prognosis and success of a complete denture. We will identify typical landmarks that we should be able to find in all patients.

During a clinical examination, anatomical landmarks that are present in nearly every patient need to be evaluated to determine if there is any distortion, abnormality, or missing landmarks due to severe alveolar bone resorption, disease processes, previous surgical alterations, or natural physical variation that would indicate a problem area.

In the maxillary arch, we should be able to identify the incisive papilla, labial and buccal vestibules, rugae, residual ridge, maxillary tuberosity, hamular notch, palatine fovea, buccal and labial frenula, midpalatine suture, and glandular area. We should then be able to determine the vibrating line that is so important to maxillary complete denture retention.

In the mandibular arch, we should be able to identify the tongue, pterygomandibular raphe, residual ridge, buccal and labial vestibules, buccal, lingual, and labial frenula, buccal shelf, retromolar pad, retromylohyoid fossa, alveololingual sulcus, lingual tubercle, and submaxillary caruncles.

Once these landmarks have been identified, we must assess how their size, shape, location, presence or absence may affect treatment and prognosis. Patients must be informed of any situation that may affect their ability to comfortably wear their complete dentures.

The denture-bearing area is that part of the attached and unattached mucosa of the edentulous ridges upon which the dentures will rest. This area becomes progressively smaller as residual ridges resorb. Maximal biting forces in

patients with complete dentures are about 5–6 times less than patients with natural or restored teeth.

Some conditions will require preprosthetic surgery and a healing period prior to the final impression being made for complete dentures.

Many complete dentures are made in mandibles with impacted wisdom teeth and the patient may experience no trouble; however, the patient needs to be informed about the risks involved with leaving an impacted tooth in place. As the ridges resorb, the bone overlying a tooth remnant or impacted tooth will resorb away, and this area will eventually only be protected from the denture by a thin layer of mucosa. When this happens, an infection or traumatic ulcer may develop. Retained mandibular third molars commonly dehisce. Have them removed!

A torus is a benign outgrowth of bony tissue covered by a thin layer of mucosa. Some maxillary tori can be left as they are, and a denture placed over them if they are not too large and do not adversely affect the retention or function of the maxillary complete denture. This may be possible because a maxillary complete denture has a broad denture-bearing area for support in the palate. Often this is managed by placing a palatal relief chamber over the torus as papillary hyperplasia will not form over these tissues, and if these areas become traumatized, they will be slow to heal due to the limited vascularity in that area.

Mandibular tori are sometimes left in place for a tooth-borne removable partial denture (RPD). When the patient loses posterior teeth and a distal extension of the flange of the denture becomes necessary, they may fail to see the need for surgical reduction, but if these are left in place, the denture base will erode the overlying soft tissues severely, causing intense pain. A tissue conditioner may help temporarily, but it will not be a long-term fix. When immediate insertion dentures are placed over tori, patients often must suffer until the tissues heal. A second surgery may not be

recommended while the tissue is thus inflamed or not completely intact.

Exostoses are bony outgrowths on the alveolar ridge. They also pose a considerable problem for a complete denture patient. The tissue is stretched tightly over them and the undercuts they have would prevent a peripheral seal for a denture from being evenly remotely possible. Food debris generally will collect in the overhangs they provide, and denture movements will denude the soft tissue overlying them.

Tori and exostoses often show signs of trauma because of their thin epithelium and prominent profile. This is particularly true after eating a hard-crusting food such as pizza. Healing will take place slowly in these areas because of the lack of adequate vascularization. Infection will have easy access to the underlying tissues.

Inflammatory fibrous hyperplasia begins as a traumatic ulcer secondary to an ill-fitting denture flange and develops into a callous-like fibroma called an epulis fissuratum. These need to be removed and the existing denture relined with tissue conditioner to promote healing. Healing must be completed before making a final impression for a complete denture.

Inflammatory papillary hyperplasia can also occur under a denture. It can arise with or without the presence of *Candida albicans* and is caused by an ill-fitting denture, wearing the denture at night, and/or poor oral hygiene. It appears as flattened or grape-like clusters, dependent upon the pressure of a denture over the area. This often occurs under a palatal relief chamber that was placed to increase suction under a maxillary complete denture.

As resorption progresses, the maxilla shrinks upward and inward, while the mandible shrinks downward and outward, leading to more and more of a posterior crossbite. As the bone resorbs, the area involved becomes less able to tolerate the presence of a denture overlying it, due to the decreased surface area and a resulting increased instability. In a severely resorbed mandible, the inferior alveolar nerve may lie on top of the residual ridge. Any pressure on this area will be painful.

A sharp mylohyoid ridge will press against the overlying tissue and make it very tender to any pressure of even a well-fitting denture. Trauma to this area can be minimized or prevented altogether if the patient inserts the mandibular complete denture in the posterior first, then slides the denture forward and down over the anterior ridge. This is an important area, as it provides an undercut that will improve retention of the mandibular complete denture.

A flabby residual ridge is simply soft tissue that becomes soft and flabby over the knife-edged bone of the alveolar residual ridges when an advanced residual ridge

resorption occurs. It is not a good support area for a denture. This condition most often develops in the maxillary anterior region because the tongue protects the lower anterior teeth from decay. A mandibular RPD is much more stable than a mandibular complete denture (CD), so the lower anterior teeth are retained to provide retention for a mandibular RPD. This results in Kelly syndrome, otherwise known as combination syndrome.

Combination syndrome or Kelly syndrome defines a situation that develops when a steep anterior guidance angle allows minimal or no posterior contacts in working, balancing or protrusive relationships when the patient goes through the chewing cycle. This results from the denture tipping anteriorly during the chewing cycle, compressing the mucoperiosteum of the premaxilla, leading to resorption of the bone of the premaxillary area. The negative pressure in the posterior can lead to a maxillary tuberosity protuberance.

Resorption can be so severe as to require augmentation with bone grafts in order to prevent idiopathic fracture of the mandible. There are cases on record where a patient has suffered a fracture from simply falling asleep with their hand on their chin. A mandibular CD is very seldom made opposing the maxillary restored arch due to the almost impossible task of achieving bilateral balance of the denture teeth against the restored teeth and the increased forces that a patient can generate while chewing against the natural or restored teeth. This can lead to rapid resorption of the mandibular alveolar ridge.

The retromolar pad is a cushioned mass of tissue, frequently pear-shaped, located on the alveolar process of the mandible behind the area of the last natural molar tooth. It is composed of nonkeratinized loose alveolar tissue covering glandular tissue, fibers of the buccinator muscle, the superior constrictor muscles, and the pterygomandibular raphe, and the terminal part of the tendon of the temporalis muscle. Since there were no teeth beneath this area, it is usually the most stable area of the mandibular edentulous ridge and should be covered by the denture flange. Trauma from ill-fitting mandibular complete dentures or the patient failing to wear a lower complete denture when chewing can also lead to this area being soft and loose. The alveolar process, on the other hand, has had teeth which have been removed and although these areas have filled in with reparative bone, it is not as resistant to resorption.

Some severely resorbed ridges will manifest as knife-edged bone under a ridge of soft tissue that feels firm. The bone under these ridges is so small that it resembles a knife in a sheath. It is very susceptible to spontaneous fracture

and will require augmentation with a synthetic material, cadaver bone or autogenous bone graft.

House classified the posterior palate by the shape of the soft palate according to how it drapes down posteriorly, relevant to developing the posterior palatal seal in a maxillary complete denture. In his classification system, a patient with a Class I palate will be able to tolerate a complete denture easily because this palate offers the broadest range of area (5–10mm) in which to place the posterior palatal seal, but it is also the hardest to locate the exact area because of its flat curvature. A Class II posterior palate is the most common in the Caucasian population, with a range of 3–5 mm. It falls between the limitations of the Class I and the Class III palate. A Class III posterior palate has a vibrating line that is easiest to locate as it drops down suddenly, but it offers only 1–3 mm of area to place the posterior palatal seal. A patient

with a Class III palate will have the hardest time of the three in tolerating a maxillary complete denture.

Neil classified the lateral throat form by measuring the height of the lingual vestibule in the retromylohyoid region. A patient with a Neil Class I lateral throat form will have over ½ in. of depth provided and is very favorable for mandibular denture retention and stability. A Neil Class II lateral throat form falls between Class I and Class III and is less than ½ in. in depth. A Class III lateral throat form has no vestibular depth and would have an unfavorable prognosis. Generally, around three-quarters of patients will have Neil Class I lateral throat form, about one-fifth will be Class II, and only 5–6% will be Class III. The lateral throat form is bounded anteriorly by the mylohyoid muscle, laterally by the pear-shaped pad, posterolaterally by the superior constrictor, posteromedially by the palatoglossus, and medially by the tongue.



## 2

## Preliminary Impressions

A preliminary impression is a negative likeness made for the purpose of diagnosis, treatment planning, or the fabrication of a tray. It is made in a stock tray that is selected from an assortment. These trays are made to fit the majority of patients, but really do not fit anyone well. They simply provide a vehicle to carry the impression material to the patient's mouth for that first impression as a starting point. This is the first impression of the patient's mouth that is made for diagnosis, treatment planning, and then for making the custom trays.

There are primarily two types of materials used for preliminary impressions: modeling compound and hydrocolloid. Modeling compound is used almost exclusively for edentulous patients. It is a thermoplastic material that is color coded for working temperature. *Thermoplastic* means that it is softened when heat is added and becomes firm when it cools. It comes in cake and stick forms, with the stick form most commonly used for border-molding a tray.

Green compound has a working temperature above 123 °F and red compound has a working temperature above 132 °F. Compound must be tempered in a water bath to obtain the correct temperature and not burn the patient. The hot water tempering bath is set at 140 °F for green compound to allow for enough plasticity to be transferred from the bath to the mouth and the compound be formed to the desired shape.

The second material used for preliminary impressions is an irreversible hydrocolloid (alginate). It is a hydrophilic gel made from seaweed, calcium sulfate, and water. *Hydrophilic* means that it mixes completely with water and it provides a very good impression in the moisture of the mouth. The setting time for alginate impressions is varied by the addition of sodium phosphate (a retarder). Fillers of diatomaceous earth and silicate powders determine the hardness of the alginate.

The dentist must hold the tray still the entire time until the alginate is set or distortion will result. A patient should never be asked to hold the tray themselves. The temperature and amount of water should be measured accurately,

or the efficacy of the material will be greatly reduced. Room temperature water (70 °F) should be used for best results.

There are three types of impression trays used for containing the impression material: Styrofoam, plastic, and metal. A Styrofoam tray is not rigid, is bulky, and does not retain the impression material well. Some plastic trays are disposable and are not very rigid. They require adhesive to retain the impression material. Reusable plastic trays are rigid and retain the material by retention holes and rim locks. Metal trays are rigid. Some are perforated and retain the impression material with retention holes. Others are solid and very rigid and retain impression material with rim lock retention. Some metal trays can be minimally adapted by bending if the tray is not made of stainless steel. Stainless trays are not adaptable.

Regardless of the tray selected, a suitable tray for a preliminary impression must conform to the general shape and size of the ridge. The flanges of the tray are contoured to follow the depth of the peripheral roll, best with a periphery wax made specifically for that purpose, although red rope wax is often substituted. The tray should allow for 3–5 mm of impression material thickness between tray and edentulous tissues. If a tray is not easily adaptable, use a smaller tray and use compound to extend and modify the tray.

With metal trays that can be adapted, use pliers or bend them with fingers. Shorten the border areas so that adaptation can start with the palate. Smooth and polish the modified areas with the pumice-impregnated rubber wheel.

A preliminary impression should have the following characteristics. It should show good borders and completely fill the peripheral roll. It should be centered in the tray and be of uniform thickness. It should show good tissue adaptation and be free of bubbles and voids. It should extend posteriorly to record the pterygomaxillary notch and the vibrating area in the maxilla. It should extend posteriorly to record the retromolar pads and the lateral throat form of the mandible.

Impression materials for preliminary impressions should be economically priced. Prepackaged impression materials are always more expensive than bulk packaged materials. Bulk packaged materials should be kept in air-tight and moisture-tight containers in a cool area. They should be kept in a relatively small container so that the material can be used in a fairly short period of time.

When a preliminary impression is poured, a double-pour technique is used. An alginate impression should be separated from the cast about 45 minutes after pouring to prevent damage to the cast as leaving stone in alginate overnight will result in extreme shrinkage and binding of the impression material to stone models. When compound is used for the preliminary impression, remove all compound from the cast before proceeding. Use a boil out tank or a hot water bath to remove wax from the cast. You must temper the cast for five minutes in warm water before a boil out tank is used, or the hot water will fracture the cold cast.

Before using the model trimmer, trim the land area of the cast with the laboratory knife. Never use a red handle knife or a scalpel for this as the blade can break and severely cut your hand. Remove all foreign material (alginate, compound, impression material, and wax) before using the model trimmer. Always have the water turned on while operating the model trimmer. Light pressure should be used to push the cast against the model trimmer wheel, so you do not lug the motor down. Slowly move the cast back and forth to achieve maximum efficiency. Brace your fingers against the table and keep hand instruments away from the model trimmer to prevent damage from the rotating wheel. Moisten the cast under running water prior to grinding on the trimmer to prevent the gypsum/water spray from adhering to the cast like concrete.

Preliminary casts should be trimmed to certain dimensions. Trim the edentulous cast to an outline form that is flat in the posterior and arched on the facial. Always trim the surface of the land areas first to visualize the width of the land area. Establish the 2–3 mm depth to the land area. (You may find it helpful to use a pencil to outline the width of the land areas before going to the model trimmer.) Trim the base of the cast to make the residual ridge parallel to the base and the bench top. Also make the cast 12–13 mm thick to a flat plain occlusion table and a flat bench top. Make the posterior side of the cast flat. Finalize the sides of the cast and the land area width. Rinse the slurry water and debris off the cast immediately. Do not use a heavy scrub brush to clean the cast or it may be damaged.

All edentulous casts must meet the same criteria. The median foundation plane is horizontal and parallel to the base of the cast when the cast is on the bench top. The cast is 12–13 mm thick in its thinnest portion. The land area is 1 mm above the tissue surface between the pterygomaxil-

lary notches and 3–4 mm in facial and 4–5 mm in posterior sections. The peripheral roll is 2–3 mm in depth from the land area. The facial outline form of the land area should follow the outline form of the peripheral roll. The sides of the cast are perpendicular to the base. The cast surface is free of voids and no nodules of stone are present. The cast is free of foreign material and debris.

## 2.1 Making the Preliminary Impression

Preliminary impressions are made in stock trays which are not capable of making an accurate impression of the vestibules. The vestibules will always be overextended on an impression made in a stock tray. The best you can do is to pick a tray that will fit over the arch and confine the impression material to that area.

An appropriate tray is selected to fit the patient as accurately as possible. Remember that a stock tray is not meant to be a perfect fit. A main concern is that the tray must be able to be placed in the mouth while loaded with alginate and removed easily as a unit without damaging the impression or the patient. The trays are made specifically for edentulous mouths. A tray made for patients with teeth should not be used to make impressions of an edentulous mouth.

If the patient has an existing denture, it can be used to determine the best size tray to be used for the impression. It is placed in an edentulous tray to determine the appropriate size tray to use for the impression (Figure 2.1). There should be adequate clearance between the denture and the tray to allow the tray to fully seat in the patient's mouth without impinging upon the tissue surfaces when it is loaded. If the patient does not have an existing denture, the dentist can spread the index finger and the second finger apart and use them to measure the width of the arches then place them in the tray to determine the best fit.

It is very important that the tray will support the alginate in the two most critical areas: the retromolar pad area of the mandible and the hamular notch (pterygomaxillary notch) of the maxillary. These areas must be captured in the impression because they are the most stable areas of the residual ridges. The tray should correspond to denture-bearing surfaces. Select a tray that is at least 3–4 mm larger than the residual ridge. The retromylohyoid space must also be captured, because these areas have such a profound effect on the retention and stability of the mandibular denture by providing a distal undercut that will resist the denture dislodging posteriorly when the patient bites on something anteriorly.

Flange areas and the palate can be modified for a more comfortable and uniform fit. Periphery wax is the best




**Figure 2.1** Place a denture in an edentulous tray to determine the appropriate size tray to use for the impression.

material used to modify and support the tray in these critical areas, but the main objective of placing the tray in the mouth and removing it without distorting the impression must not be compromised. This blue periphery wax is a mouth-temperature wax and will mold to fit the anatomy of the mouth. Red rope wax is often used instead of periphery wax but is a poor substitute as it does not have the same properties as utility wax. The impression tray is more comfortable for the patient when a strip of periphery wax is placed around the entire flange of the tray.


Border-molded trays should be comfortable and accurately cover the extensions of the arch. Blue periphery wax is applied to the complete circumference of the border of the tray (Figures 2.2 and 2.3). This wax can be modified to fit the mouth more accurately because it will become softer at mouth temperature and will allow a certain amount of border molding. The frenulum areas require special attention so that the metal of the tray does not show through.

After the periphery wax is placed on the tray, it is warmed in a warm water bath (Figure 2.4). This will soften it so that it easily molds to record the vestibular extensions in the patient's mouth. The warmed periphery wax is muscle-trimmed in the mouth to border mold the tray. When the border molding is complete, the tray is soaked in cool water to set (harden) the wax. A tray must be cool when loaded with alginate impression material before it is inserted in the mouth because the alginate will set quicker in a warm tray, which will prevent it from making a good impression.

A mouth mirror is used to pull the cheek away from the tray (Figure 2.5) to allow it to be eased into the mouth. This can sometimes be very difficult to do because some

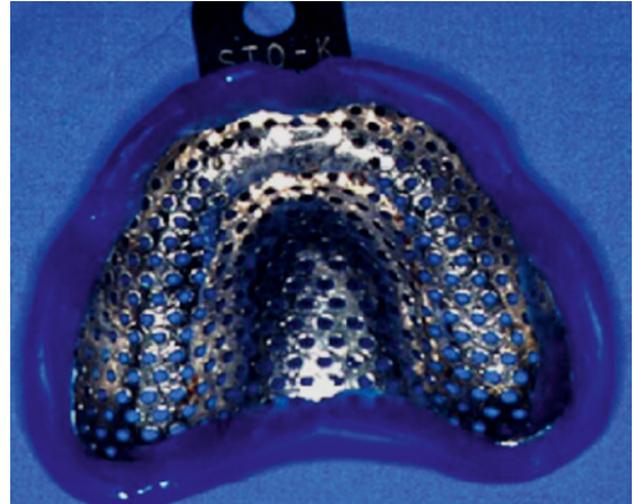

**Figures 2.2 and 2.3** Apply blue periphery wax to the tray border for patient comfort.

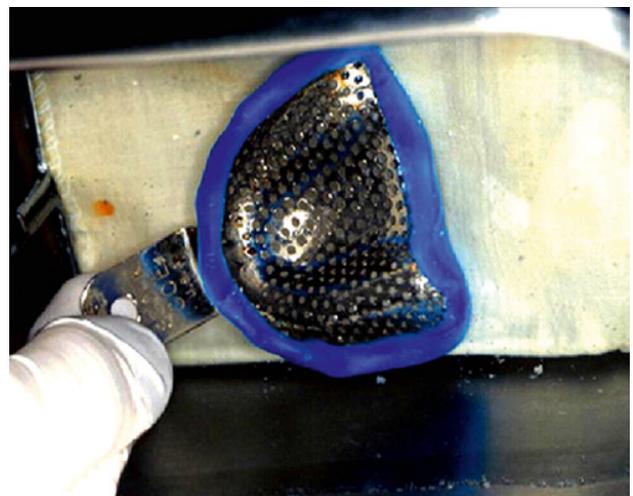

**Figure 2.4** Warm the tray in a warm water bath.



**Figure 2.5** Use a mouth mirror to pull the cheek away from the tray.



**Figure 2.7** Coat a stock tray with an alginate adhesive aerosol spray.



**Figure 2.6** Materials used to make a preliminary alginate impression.

patients' lips may have become stiff with old age. The side of the tray that is away from the mirror is inserted first, and used to push the cheek over to that side of the mouth while the fingers or a mouth mirror pull the cheek away from the other side of the tray as that side is rotated in.

Sodium alginate is by far the most popular material used to make preliminary impressions. Alginate impression material is also called irreversible hydrocolloid because it is a one-time use item. It is a very suitable impression material for this purpose, because it has enough body when it is placed in the mouth to not run excessively from the tray and sets rather quickly. It is also capable of making accurate impressions in various thicknesses of material in ill-fitting trays. Alginate impressions will pull from undercuts with a lower risk of tearing and have a very good elastic

memory to return to the arch form relatively quickly. Alginate is inexpensive when compared to most other impression materials. Major drawbacks are that the impression must be poured within 30 minutes to be accurate and the impression tray must be supported by its handle while the stone is setting.

Figure 2.6 shows a typical set-up of the materials used to make a preliminary alginate impression.

An alginate tray adhesive is applied to the inside of the border-molded tray. A thin layer of adhesive is sprayed or painted over the tissue surface of a tray. An aerosol spray can be used to coat the stock tray with alginate adhesive (Figure 2.7). The key to applying the spray adhesive is to use a thin coat. If too thick a coat is applied, it will take too long to get tacky. If alginate is placed in a tray that has wet adhesive, the alginate impression will separate from the tray when the impression is removed from the mouth. An impression will be distorted if it separates from the tray and is pressed back into place.

The correct amount of room temperature (70°F) water is poured into a vial supplied by the manufacturer of the alginate and filled to the designated line marked on the vial. The measured water is poured into a rubber mixing bowl. The powder is poured into the water for the best results and stirred to incorporate it into a homogenous mixture using a round-edged flexible spatula in a vigorous manner (Figure 2.8). Sweeping strokes are used to press the mixture against the walls of the mixing bowl to force out any air bubbles and form a thicker mix. The mix will appear dry at first but do not add additional water. After several strokes, the mix will become smooth and creamy (Figure 2.9). It is recommended to use slightly less water to achieve a thicker mix when making the maxillary impression on a gagging



**Figure 2.8** Use a round-edged flexible spatula to vigorously mix gypsum powder with water to incorporate it into a homogenous mixture.




**Figure 2.9** Use sweeping strokes to make a smooth and creamy mix.

patient. This is done rather than using warmer water as the warmer mix will set before the tray is inserted. Vacuum spatulas and other mechanical mixers are very useful and do the job very well but are not necessary once a person becomes skilled at mixing the material in a rubber bowl.

As you hold the mandibular tray by the handle in one hand, load the alginate from the back on one side of the tray and distribute it around the tray, pushing all air ahead of the mix. Make sure that the tray is filled completely (Figure 2.10).

Mandibular impressions are made by placing the loaded tray in the mouth while standing in front and to one side of the patient. The fingers of both hands are used to spread the lips and the tray is centered in the mouth with the

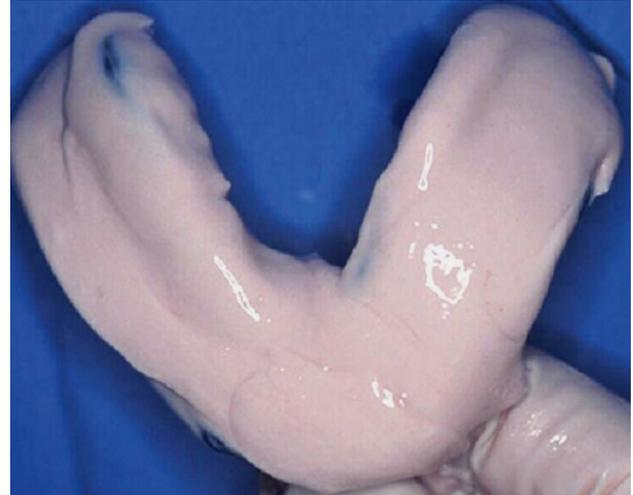

**Figure 2.10** Make sure that the impression material fills the tray completely.



**Figure 2.11** Use fingers to spread the lips and center the tray in the mouth.

handle straight out from the midline (Figure 2.11). The frenula and lips are manipulated in a procedure called muscle trimming while the tray is gently pressed down (Figure 2.12).

Have the patient raise the tongue as the tray is inserted and push it out over the handle while you press down gently, then move it to the right and left. This must be done several times while the alginate is setting, or the patient will relax and place their tongue fully back into the mouth and under the tray. The patient may worry that they will drool over themselves during this procedure, so be prepared to use a small suction tip and have paper towels



**Figure 2.12** Manipulate the frenula and lips by muscle trimming.


available close by. The lower impression should be done first because it is less likely to cause the patient to gag and they may be less anxious about making the maxillary impression.

Check the finished impression for voids and look for areas where the tray may have pressed through the alginate. This may lock the tray onto the cast and make it very difficult to remove. Be aware that alginate impressions always overextend the vestibules unless the tongue relaxes and gets under the tray in the lingual vestibule. Remake the impression if this has happened or you will not be able to make a suitable custom tray on the cast.

There should be an S-shaped curve in the lingual vestibular area on both sides of the impression. Look for the mandibular landmarks in your impression. The labial notch, labial flange, buccal notch, buccal flange, alveolar groove, and retromolar fossa should all be very identifiable in the impression. Likewise, the pterygomandibular notch, retromylohyoid eminence, lingual flange, premylohyoid eminence, lingual notch, and lingual tubercular fossa should also be easy to see (Figure 2.13).

When you “border-mold” the tray with periphery wax for a maxillary alginate impression, simply push the lips and cheeks in as you massage them (Figure 2.14). This will result in overextended vestibules which will enable you to better make a custom tray. Since the wax is very soft at this point, the tray must be removed very carefully to avoid distorting the wax (Figure 2.15). The tray is then chilled in cold water to harden the wax and slow down the setting time of the alginate.


Mix the powder and water to a thick, creamy consistency according to the manufacturer’s directions. The procedures for mixing alginate for the maxillary impression are the same as outlined for the mandibular impression. Load the

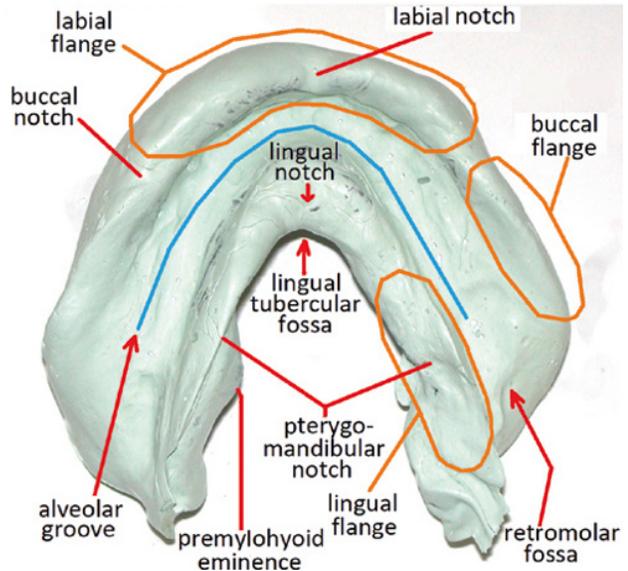

**Figure 2.13** Landmarks of a mandibular impression.



**Figure 2.14** Massage the lips and cheeks to border mold the tray.

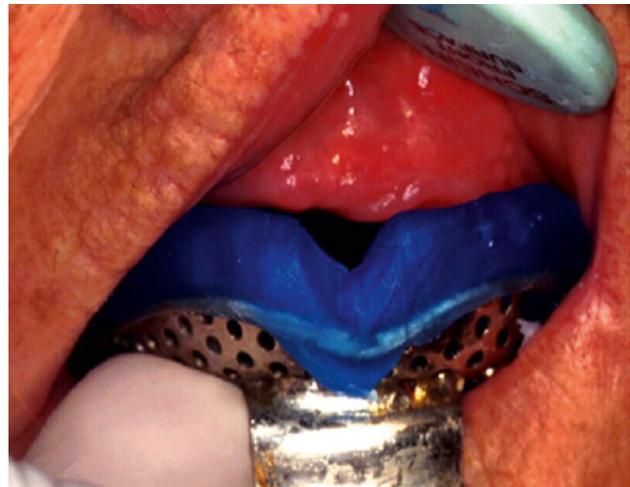

**Figure 2.15** Remove the tray carefully to avoid distorting the periphery wax.



**Figure 2.16** Load alginate mix into the tray and push air ahead of the mix.



**Figure 2.17** Use moist fingers to distribute and smooth impression material in the tray.


alginate into the tray so that the air is pushed ahead of the mix, starting at one side in the back of the tray and loading the entire tray (Figure 2.16). Distribute and smooth the material in the tray using moist fingers (Figure 2.17).

The maxillary tray is seated by the dentist from behind the patient (11 o'clock position if the right hand is used). It is placed in the mouth with one hand holding the tray handle while using the fingers of both hands to spread the lips (Figure 2.18). The tray is rotated into position, making sure that the handle comes straight out from the midline and that the tray is situated evenly over the ridge (Figure 2.19).

The anterior portion is seated first so that the air is forced out the back as the tray is seated (Figure 2.20). Any excess impression material that is forced out the back is quickly



**Figure 2.18** Seat maxillary trays while standing behind the patient.

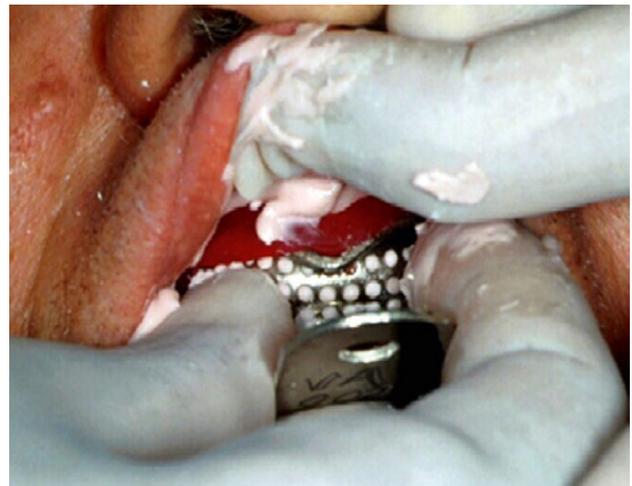

**Figure 2.19** Rotate the tray evenly into position over the ridge.



**Figure 2.20** Seat the anterior portion of the tray first to force air out the back as the tray is seated.

removed with a mouth mirror and the impression is gently border molded to form the borders of the vestibules.

If the patient does feel nauseous, have them lean forward to minimize the gagging and raise the napkin to catch any drooling (Figures 2.21). Use the mouth mirror to pull forward any excess alginate that is running out the back of the tray (Figure 2.22). Have the patient tense the stomach muscles by raising their legs while you hold the tray firmly in place for the entire time (Figure 2.23). Although the patient will feel discomfort for a short period of time, it will help to avoid having to make a second impression.

To break the seal under the denture, grasp the handle of the tray firmly, raise the cheek on one side, and remove the impression in one quick but controlled motion (Figure 2.24). You do not want to damage the patient or the impression by being too aggressive.



**Figure 2.21** Ask nauseous patients to lean forward to minimize gagging and raise the napkin to catch any drooling.



**Figure 2.22** Use a mouth mirror to pull forward any excess alginate from the back of the tray.

Check the impression for any voids that may make it necessary to remake the impression. Small bubbles (less than 5 mm in diameter) can be flicked off the stone cast but if they are larger than that, the impression should be remade. Show-through of the periphery wax does not mean that a new impression needs to be made (Figure 2.25).

Evaluate the impression to ensure that all the maxillary landmarks are captured (Figures 2.26 and 2.27). The labial flange, labial notch, alveolar groove, buccal notch, median palatal groove, buccal flange, and coronoid contour should all be very visible in the impression. Likewise, the incisive fossa, rugae, tubercular fossa, fovea palatinae,



**Figure 2.23** Have the patient raise their legs to fight nausea by tensing their stomach muscles.



**Figure 2.24** Grasp the tray handle firmly and raise the cheek on one side to break the seal and remove the impression tray in one quick, controlled motion.

## Index

### a

AccuFilm® II Double-Sided Articulating Film 105, 151, 186  
 adhesives 5, 8, 39, 45, 57, 58, 131, 199, 213, 214, 229, 245  
 air barrier coating 51, 53  
 ala-tragus line 93, 111, 126, 245  
 Al-Cote®, a tin-foil substitute 223–225  
 alginate 5–12, 205, 245, 250  
 Aluwax™ Bite and Impression Waxes 67, 68  
 alveolar groove 10, 12  
 alveolar process 2  
 Angle's classification system 246  
 anterior guidance 2, 73  
 articulating paper 108, 109, 151, 187, 191  
 articulators 59–69, 71, 73, 89, 100, 105, 106, 113, 116, 121, 130, 131, 151–153, 155, 156, 165, 169, 175, 177, 178, 182–187, 189, 197, 198, 241, 242, 245, 267  
 augmentation 2, 3, 210, 232  
 autogenous bone graft 3  
 auto-mix impression system 41  
 auto-mixing instruments 41

### b

base plates 17, 18, 21, 51–60, 62–64, 66–69, 76, 81, 82, 85, 86, 93–95, 105, 115–117, 121, 123, 126, 129–131, 140, 144, 151, 175, 234, 235, 243, 245, 246, 249, 250  
 basic tooth forms: (square, triangular, square tapering, or ovoid) 76, 215

ball abutment 229–232, 234  
 ball abutment replicas 232, 234  
 ball abutment transfers 229, 234  
 ball attachment transfers 231  
 bar system 232, 233  
 bench lathe 17  
 “bimanual technique” 130  
 bilaterally balanced occlusion 85–92, 134, 194, 199  
 bite fork 59, 60, 63  
 bite registration 59, 66–69, 71, 105, 106, 130, 175, 176, 182, 234, 245, 246, 267  
 Blu-Mousse® 67, 68, 130  
 boil out tank 6, 141  
 border molding 5, 7, 17, 20, 23, 25, 26, 28–38, 41, 44, 130, 159, 244  
 boxing an impression 47–50  
 Breaking the seal  
 “brush-bead” method 215, 218, 220, 222, 224  
 buccal flange 10, 12, 44, 172  
 buccal notch 10, 12, 171, 172  
 buccinator muscle 2  
 bulk-packaged impression material 6  
 Bunsen burner 17, 18, 25, 31, 51, 54, 55, 59, 116  
 butterfly pattern 45, 135

### c

*Candida Albicans* 2, 209  
 cap attachment transfers 229, 230, 232, 234, 235  
 Care & maintenance  
 Ceka Attachment System

centric relation, centric relation records 63, 66–71, 99, 105, 106, 108, 112, 113, 125, 129–131, 134, 151, 152, 169, 175–177, 179, 182–184, 186–191, 193, 195, 198, 199, 206, 207, 210, 234, 245, 248, 263, 264, 267  
 cheek biting 199, 206, 207, 211  
 clinical remount 105, 127, 129, 152, 175–186, 199, 200, 201  
 COE-Comfort™ Self-Cure Chairside Edentulous Tissue Conditioner 261  
 combination syndrome 2, 129, 250  
 combination tooth forms: (tapering ovoid, triangular ovoid or square ovoid) 2, 129, 250  
 complete dentures 1–3, 25, 34, 41–45, 47, 59, 71, 74, 112, 115–123, 125–127, 129, 130, 137, 139–151, 159–163, 169–173, 175, 177, 179, 184, 197–227, 229–237, 239–251, 253, 254, 263–270  
 conditioning unit 51, 53  
 containment barrier 64, 66, 69, 70, 156, 179, 181, 182  
 coronoid contour 12  
 curing oven, curing unit 17, 20, 51, 53, 244  
 cusped teeth 99, 112, 124

### d

Deflasking 139, 141  
 Denar® Slidematic facebow 59, 184  
 dental compound 28, 244  
 dental stone  
 denture adhesive (paste, powder) 123, 131, 199, 213, 214, 229

denture-bearing area 1, 249  
 denture repairs 215  
 diagnosis 5, 13, 239  
 diagnostic casts 14, 240–243  
 dimensions 6, 25, 55, 59, 63, 76, 79,  
 111, 113, 123, 125–128, 130–132, 134,  
 150–152, 195, 210, 223, 234, 241, 247,  
 249, 263, 266, 267  
 disclosing wax 169, 172–173, 204–206,  
 212, 213  
 “dished-in” appearance 111  
 disinfectant, disinfecting an impression  
 13, 250  
 distortion 1, 5, 14, 41, 51, 54, 62, 148,  
 150, 169, 257  
 double-pour technique 139  
 drying rack 14  
 DuraLay Inlay pattern resin 215

**e**

elastomeric bite registration  
 materials 67  
 Epulis fissuratum 2, 199, 209  
 equilibration 175, 187, 191, 193, 200,  
 201, 210, 245  
 esthetics, “esthetic indicators” 59, 79,  
 99, 103, 111–113, 123, 125–127, 130,  
 132, 134, 180, 204, 208, 209, 234, 241,  
 242, 247, 250, 263, 266, 269  
 exostoses 2

**f**

facebow, facebow registration 59–63,  
 65, 68, 105, 151, 152, 175, 184, 245  
 facebow transfer 59–62, 68, 105, 175,  
 176, 245  
 festooning 79, 112, 113, 115, 127, 134,  
 159, 160, 165  
 facial esthetics 126  
 final impressions 1, 2, 14, 17–23,  
 25–39, 41–45, 47–50, 204, 241, 245  
 finishing 92, 123, 152  
 flabby residual ridge 2  
 flasking 139–150, 205, 250  
 flasks 139–150, 247, 249, 250  
 flat-plane teeth 111, 112, 123,  
 124, 129  
 fovea palatinae 12–13  
 Fox occlusal plane analyzer 58  
 frenulum, frenula 1, 7, 9, 10, 25, 28,  
 29, 32, 36, 37, 44, 113, 165, 210  
 fricative sounds 56, 57, 210, 247

**g**

gagging patient 8, 12, 29, 199  
 “Golden Mean” (or “Golden  
 Proportion”) 126  
 green stick compound 25, 28, 130, 244  
 grinding 6, 108, 132, 165, 195  
 gypsum 6, 9, 47, 48, 70

**h**

Hamular notch (pterygomaxillary  
 notch) 1, 6, 13, 14, 28, 69, 135, 136  
 “Hanau’s Quint”  
 Hanau™ Wide-Vue articulator 59,  
 100, 184  
 healing collar 231  
 “high-centering” 170, 208, 219  
 high lip line 56, 74, 111, 127, 246, 247  
 horizontal condylar guide 63  
 house classification of complete denture  
 patients 3

**i**

identifying dots 78  
 immediate complete dentures  
 208–210, 219, 239–251  
 impacted teeth 1  
 implants 99, 199, 201, 202, 208, 210,  
 229–237  
 implant transfer abutments 234  
 impression materials 5–9, 11, 17, 21,  
 25, 37–39, 41–45, 106, 171, 172, 231,  
 243–245, 253  
 impressions 1, 2, 5–15, 17, 21, 23,  
 25–39, 41–45, 47–51, 106, 169, 171,  
 172, 178, 204, 205, 206, 229, 231–234,  
 236, 241, 243–245, 247, 249–251, 253,  
 257, 261, 263–269  
 incisal guide table 62, 63, 150, 151,  
 177, 179, 191  
 incisive fossa, incisive papilla 1, 12, 81,  
 82, 171, 172  
 indelible transfer stick 25, 26, 30, 123,  
 173, 205, 211  
 indexing a cast 130  
 inflammatory fibrous hyperplasia 2  
 inflammatory papillary hyperplasia 2  
 insertion appointment 127, 151, 152,  
 199, 202, 243, 250  
 instability 2, 73, 76, 111, 124, 129, 130,  
 170, 202, 204  
 intaglio surface of the denture 147,  
 170, 171, 204, 206, 250

interocclusal registration, interocclusal  
 record 62, 68, 69, 152, 170, 245  
 interpupillary line 57, 58, 111, 113,  
 125–127, 245, 246, 263  
 irreversible hydrocolloid 5, 8

**j**

jaw relationship records (JRR) 182, 185

**k**

Kelly syndrome (Plummer-Vinson or  
 Paterson-Kelly syndrome) 2, 111,  
 129, 199, 209  
 “key to occlusion” 85, 88  
 knife, buffalo, green handled 70, 123,  
 135, 136, 149, 159, 169, 175

**l**

labial flange 10, 12, 36, 44, 125, 127,  
 128, 131, 172, 239  
 labial notch 10, 12, 36, 171–173  
 laboratory handpiece 17  
 lab remount 105, 151, 154, 191, 248  
 land area 6, 14, 15, 17, 18, 21, 22, 47,  
 50–53, 56, 68, 76, 81, 82, 139, 142, 175,  
 177, 182, 200, 201  
 landmarks 1, 10, 12–14, 111, 123, 127  
 lateral condylar guide 63, 64  
 lateral throat form 3, 5  
 light-cured resins 244  
 light-curing machines 233  
 lingualized occlusion 99–109, 182,  
 183, 187–190  
 lingual tubercular fossa 10  
 Lucia jig 245  
 Lynal® tissue conditioner and temporary  
 reliner 253, 256, 261

**m**

Macroglossia 201, 207  
 master casts 17, 18, 47, 50, 51, 53, 64,  
 68, 76, 130, 139, 150, 236, 245, 250  
 maximal biting forces 1  
 median palatal groove 12  
 midline fractures 199, 208, 219,  
 220, 226  
 mini-implant system 235  
 modeling compound 5  
 Model Release Agent 51, 52  
 mold chart 74–76, 123  
 monoplan occlusion 93–97, 123,  
 197, 246

mounting casts 61–66, 69–71, 138, 150, 151, 175, 177–179, 182–186, 191, 223, 246  
 muscle-trimming 9, 10, 44, 59  
 mylohyoid ridge 2, 171

## n

Neil's classification of lateral throat form 3  
 neuromuscular control 199, 209, 210, 213  
 “neutral zone” 112, 132

## o

occlusal equilibration 175, 191, 193, 200, 201  
 occlusion rims 51–56, 66–69, 81, 83, 85, 87, 111, 200, 201, 234, 245  
 “orange peel appearance” 113, 134  
 organizational chart of representative molds 76  
 overdentures, advantages, disadvantages 229

## p

packing 47, 115, 139–150, 187, 195, 205, 235, 250  
 palatal relief chamber 1, 2  
 palatal seal 3, 13, 14, 17, 21, 28, 30, 36–38, 45, 115, 116, 125, 134, 135–136, 170, 171, 199, 203–206, 211, 244, 254  
 palatal torus 45, 170, 203, 205, 206, 208  
 papillary hyperplasia 1, 2, 199, 206, 209, 210  
 patient remount 170, 199, 210, 248  
 periphery wax 5–7, 10, 12, 13  
 permanent soft liners (silicone elastomers) 213, 214  
 planes: horizontal, sagittal and vertical 6, 67  
 plaster/pumice mix 47–49  
 Playdough (Play-Doh) 215  
 polishing 159, 163, 165, 166, 224, 268  
 polishing lathe 159, 162  
 polishing materials 5, 223  
 polyether  
 polymethylmethacrylate 263  
 polysulfide rubber base 41  
 polyvinylsiloxane 41, 43, 179, 215, 217, 219, 220, 222, 233, 263, 265, 267

posterior palatal seal (arbitrary, functionally-generated) 3, 13, 14, 21, 28, 30, 36–38, 45, 115, 116, 125, 134–136, 170, 171, 199, 203–206, 211, 244, 254

pouring a final impression 47–50  
 pouring time 6, 8, 14  
 preliminary casts 6, 17, 25, 243, 245  
 preliminary impressions 5–15  
 premylohyoid eminence 10  
 processing error 51, 139, 150–153, 182–184, 187, 191, 248  
 pterygomaxillary seal 1, 2, 5, 6, 10, 13  
 prepackaged impression material 6  
 preprosthetic surgery 1  
 presses 112, 144  
 pressure indicator paste (PIP) 39, 169, 170–172, 175, 191, 211, 236, 250  
 processing  
 “processing error” 51, 139, 150–153, 182–184, 187, 191, 248  
 protrusive relation record 175  
 pterygomandibular raphe 1, 2, 10  
 pterygomaxillary notch 5, 6  
 pumice 5, 159, 175, 177, 219  
 “punched-out” appearance 111

## r

“rabbeting” 220, 223–225  
 rare earth magnets 229, 236  
 red rope wax 5, 7  
 “refractory patient” 240  
 Regisil® 67  
 registration index 69, 71  
 relief chamber 1, 2, 232, 233  
 remount casts 105, 106, 175–186, 191  
 remount index 151, 152, 155–157, 169, 175, 177, 179, 183  
 remounting casts 130  
 repairs 73, 74, 199, 203, 204, 207, 208, 211, 215, 217–227, 232  
 Reprosil® 50ml cartridges  
 resorption 1, 2, 14, 111, 129, 130, 132, 170, 199, 203, 204, 208, 213, 226, 227, 229, 242, 250  
 retarder 5  
 retromolar pad 1, 2, 5, 6, 13, 14, 21, 55, 56, 69, 76, 111, 123, 124, 126, 197  
 retromylohyoid eminence, retromylohyoid space 6, 10, 119, 202, 203  
 reusable plastic trays 5

reverse crowfoot pattern 105, 187  
 rugae 1, 12, 115, 165, 171, 210  
 rule of BULL 191  
 ruler, flexible plastic 81, 83

## s

sandpaper, wet/dry 70, 175, 177, 197, 198  
 seating the tray 44  
 selecting teeth 73  
 selective grinding for esthetics 132  
 selective pressure technique 28  
 semiadjustable articulator 100, 241, 242  
 separating impressions from casts 220, 222  
 separating medium 139, 141, 142, 220, 222  
 setting articulators 63, 67  
 setting teeth 73, 76, 99–109, 111, 206  
 setting time 5, 10, 41  
 S-shaped curve 10  
 stippling 112, 113, 115, 134, 165  
 stock trays 5, 6, 8, 17, 245  
 superior constrictor muscle 2  
 surgical stent 247–250  
 styrofoam tray 5

## t

Take 1™ bite registration material 67, 245  
 tempering bath 5, 28  
 tempering thermoplastic materials 5  
 temporalis muscle 2  
 thermoplastic materials 5  
 third point of reference 59, 63, 64  
 tissue conditioner relines 243, 250, 253, 255, 261  
 tongue exercises 253  
 tooth repairs 73, 74, 199, 203, 204, 207, 208, 211, 215  
 torus, tori 1, 2, 45, 135, 170, 171, 199, 203, 205, 206, 208, 250  
 transfer analog 234  
 trays 5–8, 11, 17–23, 25–39, 231, 243–245, 263–265, 267–269  
 treatment planning 5  
 Triad® TruTray™ VLC Custom Tray Material, Dentsply Triad Denture Base Material 17, 19–22, 51–53, 76, 77, 115, 140, 215, 234, 244  
 trimming casts 6, 47, 50, 139, 142

trimming flash 257  
trimming impressions 44  
troubleshooting 199–214  
Trubyte® Tooth Indicator-Dentsply  
    International 75, 76  
tubercular fossa 10, 12  
two-tray impressions, two-tray method  
    244, 245

**V**

vacuum spatulas 9  
vermillion border

vertical dimension at rest (VDR) 125, 234  
vertical dimension of occlusion (VDO)  
    79, 123, 125, 152, 210, 234, 247, 263,  
    266, 267  
vestibule, vestibular depth 1, 3, 6, 10,  
    12, 17, 18, 21–23, 25, 28, 29, 43, 50–52,  
    55, 56, 81, 82, 89, 111, 113, 124, 131,  
    172, 200, 204, 210, 243–245  
vibrate the impression 14

**W**

waxing instruments 47, 49, 123

wax try-in 79, 111–113, 123–134, 210,  
    241, 247

**Z**

zero-degree teeth 76, 92, 93–97,  
    99–102, 104, 112, 123, 125, 151, 152,  
    154, 169, 182, 197, 198, 200, 206, 207,  
    210, 246  
ZEST Anchors, Inc., LOCATOR®  
    Attachment System  
ZnOE Bite Registration Paste 67